

# Over Bed Table

Item # OB2845

## User Manual



Innovative Healthcare Solutions

Over bed table installation parts listed below (See Figure # 1)

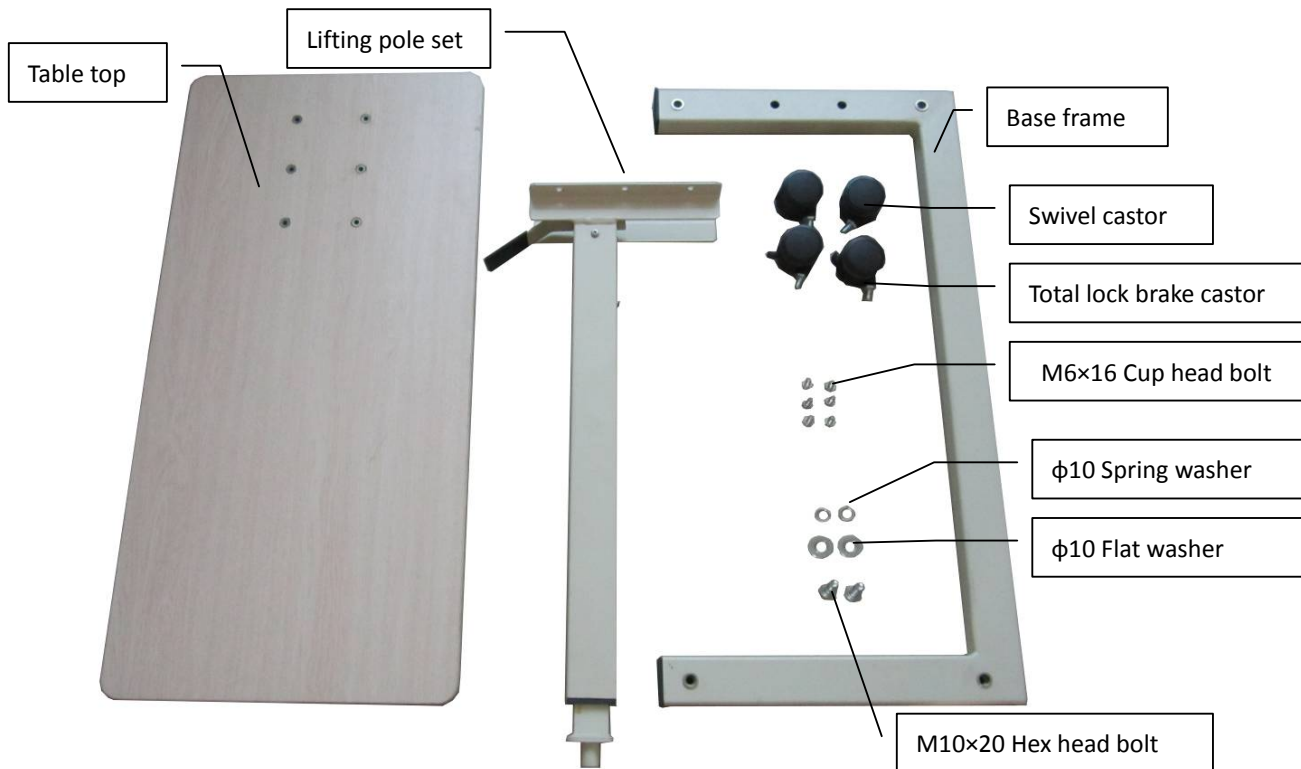


Figure # 1

### 1. Installing the table top

Place the table down face down on a flat surface with the bolt holes facing up. Match the six bolt holes on the lifting pole with bolt holes on the table top.

Fasten 6 pcs M6 X 16 cup head bolts with a screwdriver. (See Figure # 2)

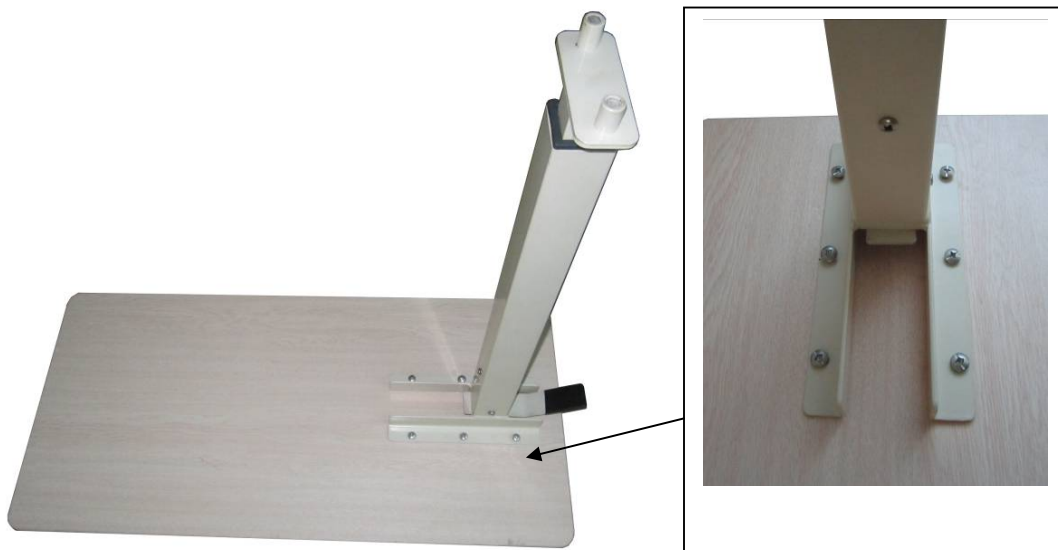


Figure # 2

## 2. Installing the castors

Install the total lock brake castors diagonally. Insert the caster stem into the M10 bolt hole in the base frame. Fasten the caster stem by using an open spanner (See Figure # 3)

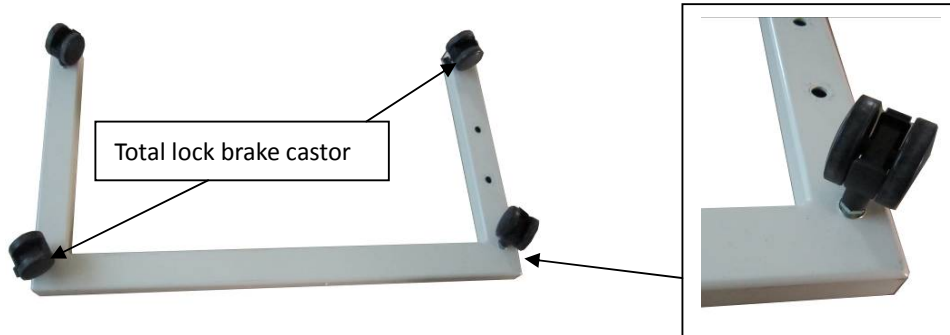


Figure # 3

## 3. Installing the base frame

Matching the  $\phi 18$  hole in the base frame with the M10 bolt hole in the lifting pole. Insert the screw nuts of the lifting pole into the base frame. Insert the  $\phi 10$  spring washer and  $\phi 10$  flat washer into each one of 2 pcs M10 X 20 hex head bolts, and fasten the bolts with a spanner. (See Figure # 4)

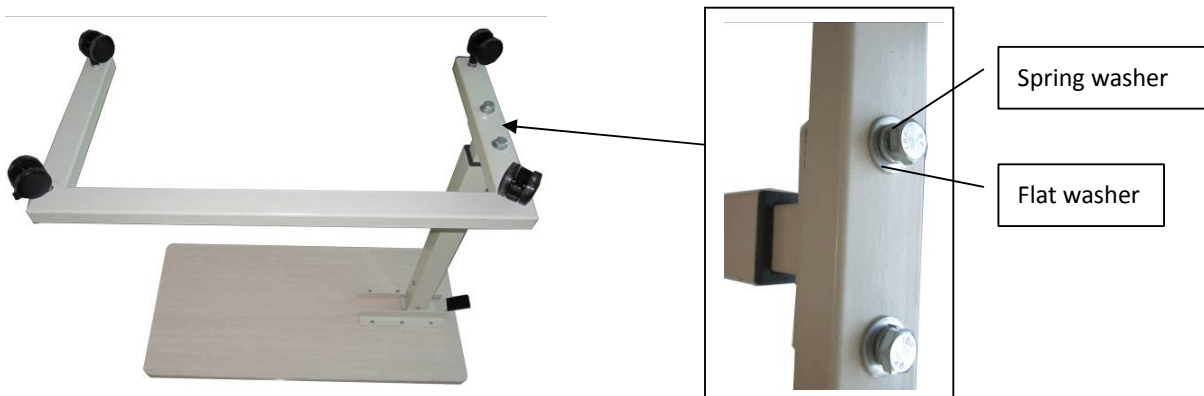


Figure # 4

#### 4. Operation of the over bed table

Place the over bed table in the upside down position and remove the self-tapping screw (for shipping purpose).

The over bed table has an adjustable height from 28" – 45"

To raise the over bed table, lift the table top with your hand (do not press the lifting handle).

To lower the over bed table, push up the lifting handle with one hand and press down the table top with the other hand. (See Figure # 5)



Figure # 5